

# REVIAN<sup>®</sup> RED SYSTEM

## The REVIAN System includes:

- 1- REVIAN Cap
- 1- Accessory Kit (AC Adaptor, USB Cable, Adjustment Headliners, Travel Bag)
- 1- Charging Support Base
- 1- Downloadable REVIAN app (mobile device not included)

## Quick Start Guide:


## Regrow in **3** easy steps!



- 1** Charge Cap with the included charger. Download REVIAN app & click NEW USER or enter LOGIN from REVIAN.com.

Follow on-boarding instructions within the app to get started.


\* Visit MANUAL in menu  for more detailed instructions.

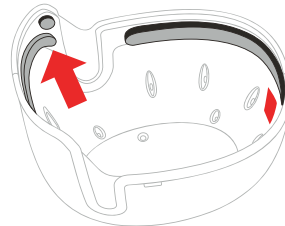



- 2** A front sensor awakens the Cap. The Cap must have a snug fit against the forehead to operate.


Remove & replace **only the back headliner** to adjust fit.

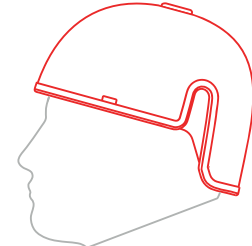\* To test fit, visit MANUAL in menu  - Adjusting the Cap



- 3** Place the Cap on your head and press the yellow  icon on the app Treatment page. Next, under Settings press **connect** to REVIAN RED.

Once connected **press**  to begin a daily 10-minute treatment.


 Cap Bluetooth Connection will only be seen in the app, not in the phone settings.

Disconnect from any current Bluetooth devices before connecting to the Cap.


## Charging steps:



To charge, **attach the USB cable** to the **AC adaptor**, then plug the adaptor into an outlet. Place the Cap on the charging base and **open the USB cover** located on top of the Cap.

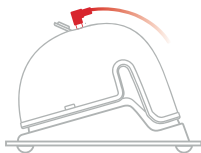

**Insert** the micro-USB **cable** into the micro-USB **port** located on top of the **Cap**. The charging port will **illuminate blue** when successfully charging.

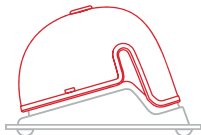

Once the Cap is **fully charged**, the charging port **illumination will turn off**. **Unplug** the micro-USB **cable** before **removing** the Cap from the support base. After treatment, place the Cap back on the support base for storage.

## IMPORTANT!

- For best results, **charge Cap 60 minutes** before use.
- The front headliner should be replaced for sanitary purposes only.
- **Disconnect** from and previously connected Bluetooth devices, such as smart watches prior to use.
- A treatment can be paused by pressing the pause button on the app or simply by removing the Cap.

For assistance, please visit us at:  
[REVIAN.com/support](https://www.revian.com/support)

1240-000005\_01

**REVIAN<sup>®</sup>**  
SCALP | SCIENCE

Connect with us at  
[REVIAN.com](https://www.revian.com)

